

**DATE:** 5/16/2011

**PRODUCT:** NeuraJet17

**TOPIC:** Log Form Creation for Print Driver

**ISSUE OR PROBLEM:** Log forms must be set up for the NeuraJet17 print driver to match actual Paper Log Size; sometimes the automated process fails and this must be done manually

**SOLUTION:**

Standard log forms are normally installed during Print Driver setup. These forms may also be setup manually.

Log 11 Inch Wide Paper	Log 14 Inch Wide Paper
Log 17 Inch Wide Paper	Log 8.75 Inch Wide Paper

**RESOLUTION/STEPS:**


After Installing the NeuraJet17 Print Driver, select the NeuraJet17 printer through your Windows Control Panel Devices and Printers.

With the printer selected, choose *Print Server Properties* (accessed either at the top of the window or through a right click menu).

On the Forms tab, select *Create a new form*.

Enter the Form name and enter values shown.

Select OK when done. These forms will now be selectable in the printer's page setup tab when you print. Select the correct form for *Page Size*: to match the width of the paper you have loaded into your NeuraJet17.

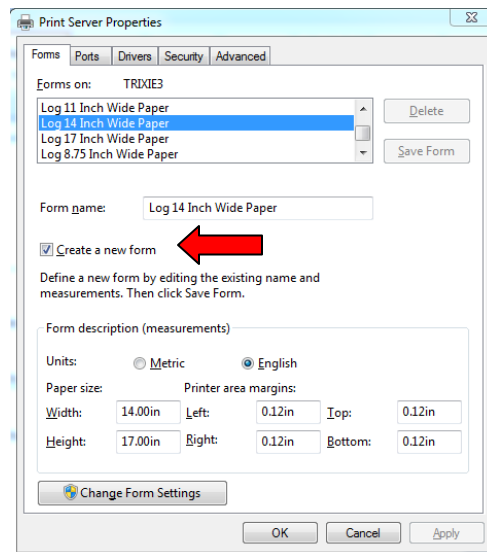

Form Name	Width	Height	Left/Right/Top/Bottom Margin
Log 11 Inch Wide Paper	11 in	17 in	0.12 in
Log 14 Inch Wide Paper	14 in	17 in	0.12 in
Log 17 Inch Wide Paper	17 in	17 in	0.12 in
Log 8.75 Inch Wide Paper	8.75 in	12.50 in	0.12 in

**NOTE:** Once these Log forms are installed and shared on a print server, any clients accessing the shared printer will automatically have access to the forms.